



# Patient Lift Sling Type Reference

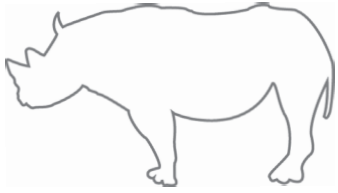


Tuffcare Rhino slings are made with a three layer soft but sturdy polyester fabric as well as a polyester mesh type of fabric which provides comfort and protects patients who may have sensitive skin. Along with preventing the patient from sliding in the sling, our Tuffcare Rhino slings are universal and will accommodate any 2, 4 or 6 point spreader bars from other suppliers.

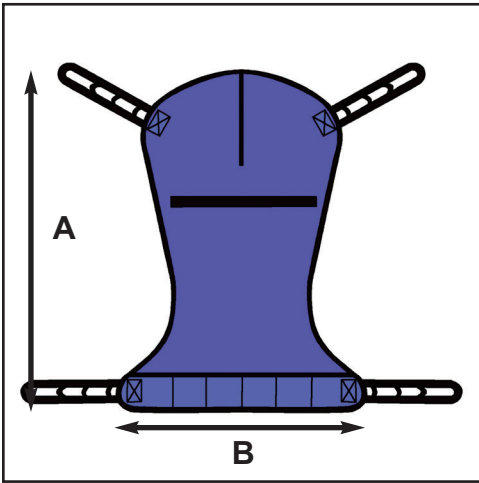

## Full Body Sling

Tuffcare full body slings offer both total head and neck support along with a padded lower section that provides additional under thigh comfort and support. These full body style slings accommodate the whole body to be supported in the sling with the arms remaining inside the sling straps and can be used for amputee patients in a recline position. The patient legs will hang on the outside of the full body slings from the knee to the foot area.

Model #	SIZE	A	B	Recommended Weight
P110	S <span style="color:blue">■</span>	47"	39"	55-110 lb.
P111	M <span style="color:green">■</span>	53"	43"	99-210 lb.
P112	L <span style="color:magenta">■</span>	58"	47"	198-350 lb.
P113	XL <span style="color:black">■</span>	63"	49"	270-600 lb.

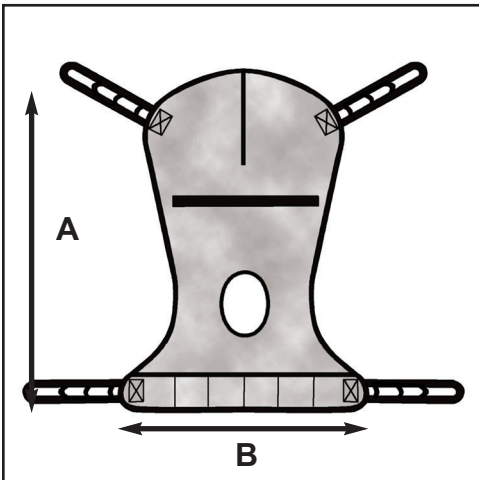


## Full Mesh Body Sling with Commode Opening

Tuffcare full mesh body slings with commode opening offer total head and neck support along with a padded lower section that provides additional under thigh comfort and support. These Tuffcare full mesh body slings with commode opening will hold the body with the arms remaining on the inside of the sling straps and can be used for amputee patients in a reclined position. The patient's legs will hang outside of the sling from the knee to foot. Our Tuffcare full mesh body slings are designed for use in the shower, as the grey mesh material will quickly dry upon being splashed. This sling's commode opening is ideal for toileting usage.

Model #	SIZE	A	B	Recommended Weight
P120	S <span style="color:blue">■</span>	47"	39"	55-110 lb.
P121	M <span style="color:green">■</span>	53"	43"	99-210 lb.
P122	L <span style="color:magenta">■</span>	58"	47"	198-350 lb.
P123	XL <span style="color:black">■</span>	63"	49"	270-600 lb.

COMMODOE OPENING 9" X 11"

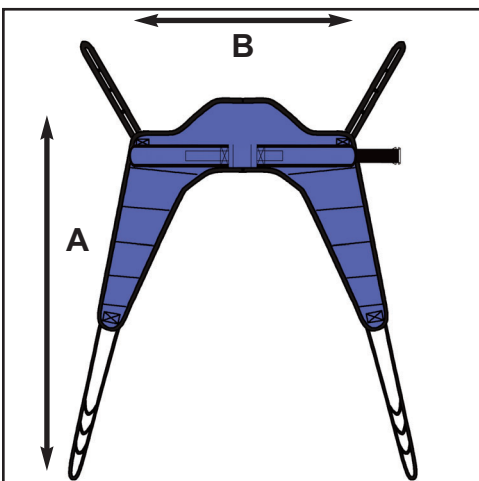


## Universal Toileting Sling

Tuffcare universal toileting slings are designed to provide easy access for toileting along with providing personal hygiene cleaning and adjustment for clothing when needed. Our universal toileting slings should only be used when you are transferring a patient who has full head and neck control while in a seated position.

Model #	SIZE	A	B	Recommended Weight
P161	M <span style="color:green">■</span>	60"	32"	99-210 lb.
P162	L <span style="color:magenta">■</span>	62"	34"	198-315 lb.
P163	XL <span style="color:black">■</span>	65"	37"	310-450 lb.


Belt: medium 22" - 40"; large 35" - 53", X large 40" - 60"


## Washing Instructions


 Machine wash at 180°F

 Low temperature dry

 Sling can suffer damage during washing and drying and should be checked carefully before each use.



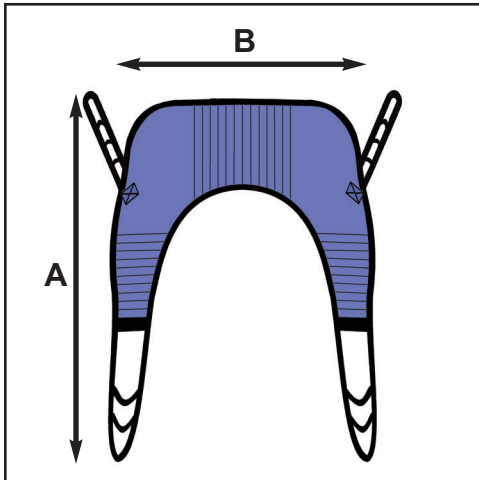

 Do not use bleach


 No dry clean



### Universal Sling



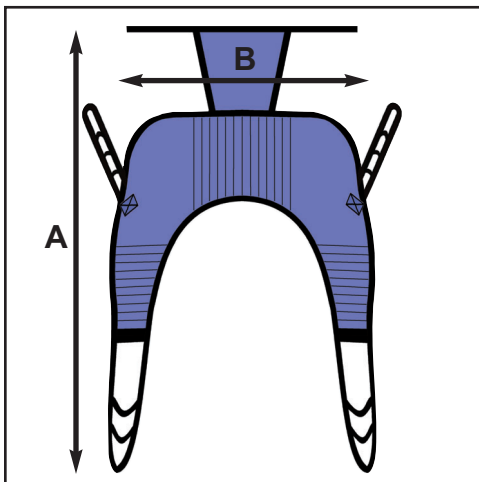

Tuffcare universal slings are designed for patients that are partially dependent and are able to hold their head up and maintain their upper body strength. The universal U sling wraps around the patient's body and this will eliminate the need to pull or drag the sling underneath the patient's body along with having a large opening to accommodate toileting when needed.

Model #	SIZE	A	B	Recommended Weight
P130	S 	62"	30"	55-110 lb.
P131	M 	65"	33"	99-210 lb.
P132	L 	68"	42"	198-350 lb.
P133	XL 	72"	46"	270-600 lb.




### Universal Sling with Headrest

Tuffcare universal slings with headrest are designed to accommodate patients that are partially or totally dependent. The universal sling with headrest has extra head support material for patients that cannot hold their head up or have very limited upper body strength. Our universal U sling with headrest wraps around the patient's body and will eliminate the need to pull or drag the sling material under the patient's body and also provides a large opening to accommodate any toileting when needed.



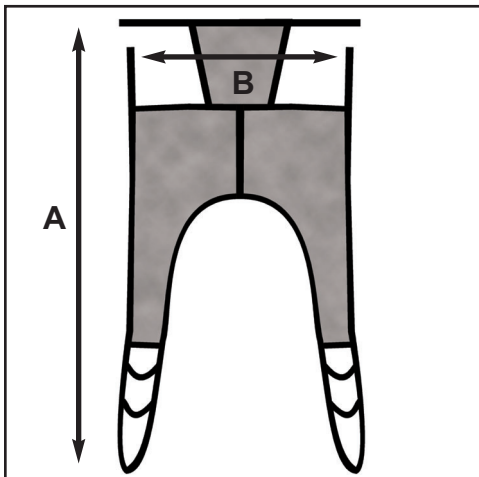

Model #	SIZE	A	B	Recommended Weight
P135	S 	75"	30"	55-110 lb.
P136	M 	78"	33"	99-210 lb.
P137	L 	81"	42"	198-350 lb.
P138	XL 	85"	46"	270-600 lb.

HEADREST: 16" / 14" W x 13" H




### Universal Sling with Headrest, for Shower

Tuffcare universal mesh slings with headrest are designed for patients that are either totally or partially dependent while using in the shower. The universal mesh slings with headrest has extra head support material for patients that cannot hold their head up or have very limited upper body strength. Tuffcare universal mesh slings with headrest wrap around the patient's body which will eliminate the need to pull or drag the sling material underneath the patient's body along with having a wide opening to accommodate any toileting when needed.

Model #	SIZE	A	B	Recommended Weight
P170	S 	75"	20"	55-110 lb.
P171	M 	78"	25"	99-210 lb.
P172	L 	81"	30"	198-350 lb.
P173	XL 	87"	36"	270-600 lb.

HEADREST: 16" / 14" W x 13" H

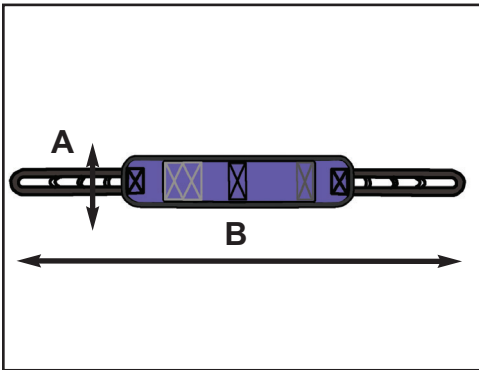


### Hip Lift Strap

The Tuffcare lower body hip lift strap works with all standard patient lifters with T type hanging bar. The hip lift strap mounts onto both of the patient's legs. The strap gently elevates the lower half of a patient's body and gives caregivers easy access for cleaning, changing, and careful delivery needs.

Model #	SIZE	A	B	Recommended Weight
P165	S <span style="color:blue">■</span>	6"	61"	55-110 lb.
P166	M <span style="color:green">■</span>	6"	64"	99-210 lb.
P167	L <span style="color:magenta">■</span>	6"	69"	198-350 lb.
P168	XL <span style="color:black">■</span>	6"	74"	270-600 lb.

Belt: small 11" -17"; medium 14" - 20"; large 19" - 25", X large 24" - 30"

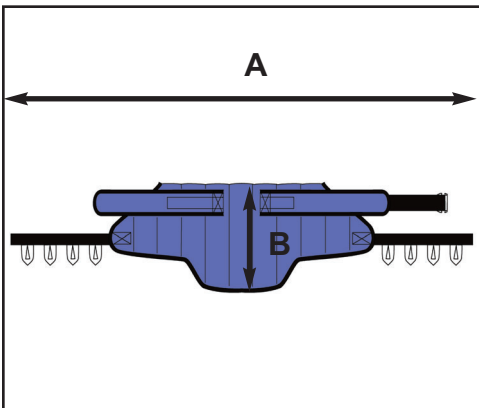


### Standing Sling

The Tuffcare standing sling is designed to be used with a stand up lift only. It will accommodate patients who have weight bearing capacity along with head and neck control, and are able to bend at the knees, hips and ankles. Our Tuffcare standing sling provides patient dependence along with effectively transferring a patient from the bed to a toilet, shower chair or standard chair.

Model #	SIZE	A	B	Recommended Weight
P141	M <span style="color:green">■</span>	68"	16"	99-210 lb.
P142	L <span style="color:magenta">■</span>	74"	16"	198-300 lb.
P143	XL <span style="color:black">■</span>	84"	16"	270-450 lb.


Belt: medium 22" - 40"; large 35" - 53", X large 40" - 60"

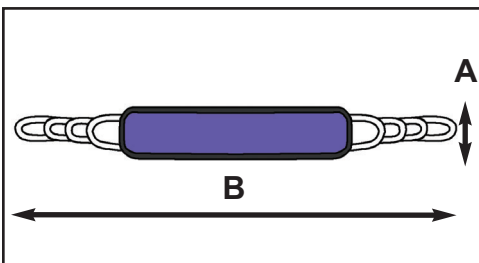

### Sit to Stand Buttock Strap

The buttock strap is to be used in conjunction with a universal sling and not to be used by itself. The buttock strap works perfect when the patient needs additional help and security while using their sit to stand lift. The buttock strap works by placing the strap underneath the patient's bottom. This will allow the patient to sit in the buttock strap like a swing, and will relieve any additional weight on their lower half.

Model #	SIZE	A	B	Recommended Weight
P147	L <span style="color:magenta">■</span>	7.5"	71"	99-450 lb.

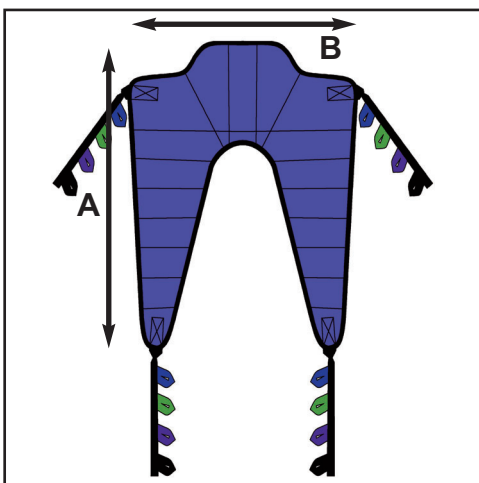

### Transport Sling

The Tuffcare Transport sling is designed to be used only with a stand up sling. It will accommodate patients who have weight bearing capacity, have full head and neck control, are able to bend at the hip, knees and ankles, and are able to sit up on the edge of the bed's surface. The Transport sling is used for simple transfers or transporting patients for longer distances while being in a seated position and will accommodate toileting due to its large opening.

Model #	SIZE	A	B	Recommended Weight
P152	L <span style="color:magenta">■</span>	40"	34"	99-230 lb.
P153	XL <span style="color:black">■</span>	44"	46"	225-450 lb.

# Patient Lift Sling Size Reference Guide

Our slings come in different sizes and have different limitations to the amount of weight that they can support. Although most slings are offered in multiple sizes, which typically support more weight as the size increases, some sizing options are more limited than others. From the chart below, based on sling type, you can correlate the user's weight and area measurement to choose the perfect sling size.

	Measure	Small	Medium	Large	X-Large
Sling Type / Model Number	Point	55 - 110 lbs	99 - 210 lbs	198 - 350 lbs	270 - 600 lbs
Full Body Slings / P110 - P123	A	≤ 39"	37" - 50"	47" - 62"	≥ 57"
Universal Slings / P130 - P138	A	≤ 37"	35" - 48"	45" - 60"	≥ 55"
Toileting Slings / P161 - P163	B	X	22" - 40"	35" - 53"	40" - 60"
Universal Shower Slings / P170 - P173	A	≤ 37"	35" - 48"	45" - 60"	≥ 55"
Hip Lift Sling / P165 - P168	C	11" - 17"	14" - 20"	19" - 25"	24" - 30"
		Small	Medium	Large	X-Large
		55 - 110 lbs	99 - 210 lbs	198 - 315 lbs	310 - 450 lbs
Standing Sling / P141 - P143	B	X	20" - 40"	35" - 53"	40" - 60"
Buttock Strap / P147	B	X	X	≤ 60"	X
Transport Sling / P152 - P153	B	X	X	≤ 53"	40" - 60"
X - No size		≤ - Less or Equal to		≥ - Greater or Equal to	

